

# Summer Internship Report

Abhishek Anil Deshmukh

July 22, 2019

## 1 Heat Shock

Heat Shock is given to cells for getting the mutant plasmids into them, We keep them in a temperature of  $-4^{\circ}\text{C}$  for 15 minutes then move them to  $42^{\circ}\text{C}$  for 35 seconds and then move back to  $-4^{\circ}\text{C}$ .

In this process all the plasmids do not enter the cells, so you incorporate a X-Resistant gene into the plasmid (Here X can be anything, we chose Ampicillin) and now after incubation of about an hour in a rich medium at 220RPM and  $37^{\circ}\text{C}$ . After that colony making is done at  $37^{\circ}\text{C}$  for atleast 12 hrs (we did it for about 18 hrs).

## 2 Ultracompetent Bacterial Cells

These are the cells which clone plasmids for large number of them. (Basically PCR for large sequences)

PCR has a limit of about 4kb, so for cloning large sequences Ultracompetent Bacterial Cells are used.

These cells are given a heat shock and then left in medium for an hour, later on they are taken for colony making for atleast 12 hrs at  $37^{\circ}\text{C}$ .

## 3 lb Broth Preparation

2.5gm of lb broth powder per 100 ml of water, the boxes of lb broth and lb agar look surprisingly similar, so watch out for that. 4mg of lb agar powder per 100ml of water. If we are making x ml of lb broth solution then the size of the container should be 5x to 6x ml as the size increases when bacterial or cells feed on it.

## 4 Competence

The ability of a cell to alter its own genetics by taking up extracellular DNA from the environment. It is brought into cell by environmental conditions like

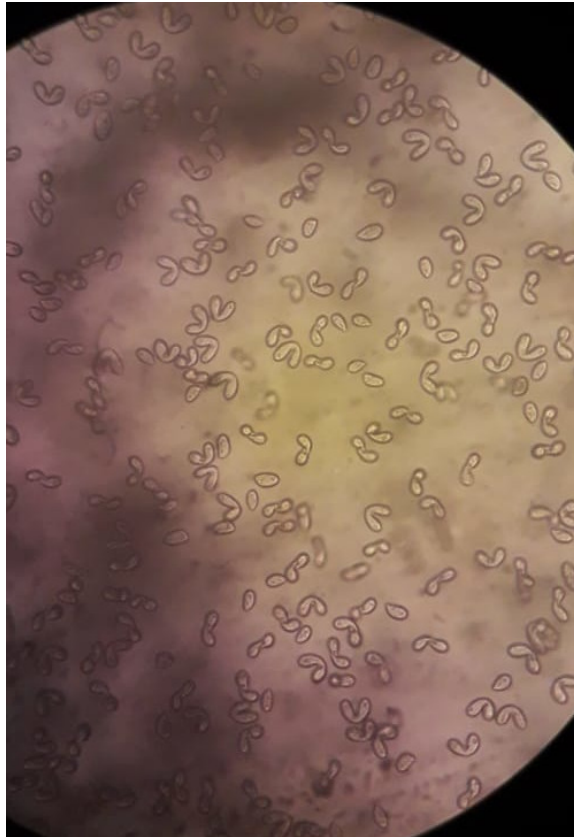


Figure 1: Conjugation

starvation. Conditions inducing sporulation often overlap with conditions inducing competence.

## 5 Transformation

It is a process in which the competent cell intakes DNA from outside and ends up changing its own DNA. Transformation usually produces a mixture of relatively few transformed cells and an abundance of non-transformed cells, a method is necessary to select for the cells that have acquired the plasmid. The plasmid therefore requires a selectable marker such that those cells without the plasmid may be killed or have their growth arrested. Antibiotic resistance is the most commonly used marker for prokaryotes.

## 6 Conjugation

At the 8<sup>th</sup> hour

## 7 Agarose-Gel

For 100ml of agarose gelcast of 0.8

1. Tape the gel caste tray
2. Put the comb on the end of the cast tray
3. Agar(0.8g)+ 1x TAE(100mL) in a flask
4. Oven it till dissolution (In this case it took about 5 minutes), check after every minute or so.
5. Let the flask cool before putting in ETB, but don't let it cool so much as to start polymerization.
6. Put in EtBr(2 $\mu$ L) shake well.
7. Put it in the tray
8. Move the bubbles away from the well, towards the ends.

## 8 Cell Fixation

Cell fixation is the process of arresting the cells completely till the protein level. This can be done in two ways.

1. PFA  
Paraformaldehyde. The sample of cells are centrifuged to precepitate, the floating part is thrown away, PFA is added quickly. After sometime the cells are suspended again, to the original concentration.
2. Cooling  
This process involves slow cooling down of the cells.

## 9 Cell Counting

Hemocytometer is a counting-chamber device orginally designed and usually used for counting blood cells.

1. Take a batch of starved cells(basically anycells which you want to find the cell count).
2. Fix them with PFA(0.5 $\mu$ L)
3. Add the solution(10 $\mu$ L/ chamber) in both the chambers.(keep the pipet horizontal for faster spreading)

4. Allow the chamber to fill(Thanks to capillary action, interference is not required).
5. View under the microscope.
6. Determine the no.of cells  
 Count the number of cells in each square, 4 squares in each chamber,count both the chambers. Divide the number by 80.(divide by 8 so we get average number per square and then by 10 because of the constant) The number which we get is that many million cells/mL.

## 10 DAPI Staining

4,6-diamidino-2-phenylindole(DAPI) staining is the process of arresting the cells and staining them with a UV Sensitive dye(DAPI).

1. Centrifuge 1ml of cell under 1.1G for 2 minutes and discard the supernatant(the floating stuff).
2. PFA 500 $\mu$ L(At the time when you want to arrest them)(In our case at 8 hours for when they were kept together).
3. Tap and Invert.
4. Centrifuge under 1.1G for 2 minutes and discard the supernatant(the floating stuff).
5. Add 200 $\mu$ L of 10mM Hepes(A buffer), tap and Invert.
6. Centrifuge under 1.1G for 2 minutes and discard the supernatant(the floating stuff).
7. Add 200 $\mu$ L of DAPI.(Preferably in dark)
8. Leave in dark for 10 minutes.
9. Centrifuge under 1.1G for 2 minutes and discard the supernatant(the floating stuff).
10. Take 5 $\mu$ l of cells and prepare slide.


DAPI is light sensitive so stored in darkness in fridge.

## 11 DMC Preparation


Be very precautious of contamination while making it.

1. Use autoclaved MQ
2. discard old DMC from the DMC deducted bottle.(each and every drop)
3. Preparation Requirements

(a)

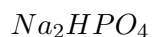

(b)

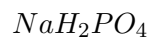

(c) Sodium Citrate

(d)




(e) MQ water


(f) renin pipette (1ml)

Procedure for Preparation:-

1. Add 1ml of

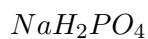

. Then 1ml

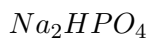

2. Add 100mm Sodium Citrate (1.7ml).
3. Add 100mm Calcium Chloride (1.5ml).
4. In measuring cylinder pour exactly 1L autoclaved MQ.
5. Remove 5.2mL of Autoclaved MQ by renin pipette and pour in DMC bottle.
6. Place cap and mix properly.
7. Autoclave it for 40 mins along with
$$CaCl_2$$
and
$$MgCl_2$$
.
8. After 12-14 Hrs it is ready for starvation.

## 12 SDS-PAGE casting

Sodium Dodecyl Sulfate - PolyAcramide Gel Electrophoresis This is gel-casting for running protein.

1. Set the caster and put water into it to test for leakages, leave for 15 mins.
2. When you are sure the caster does not cause any leaks pour out the water and get rid of excess of it with the help of tissue.
3. Make the resolving gel according to the concentration required.(5ml for every cast)

Table 1: Resolution Gel

Components	7%	8%	9%	10%	12.5%	15%
30% Acrylamide	1.16ml	1.33ml	1.5ml	1.66ml	2.08ml	2.5ml
MQ Water	2.48ml	2.32ml	2.15ml	1.98ml	1.57ml	1.5ml
1.5M Tris (pH 8.8)	1.25ml	1.25ml	1.25ml	1.25ml	1.25ml	1.25ml
10% SDS	50µl	50µl	50µl	50µl	50µl	50µl
10% APS	50µl	50µl	50µl	50µl	50µl	50µl
TEMED	5µl	5µl	5µl	5µl	5µl	5µl

Table 2: Resolution Gel(5ml)

Components	1 Stack (2ml)%	2 Stacks (4ml)	3 Stacks (6ml)	4 Stacks (8ml)
30% Acrylamide	340µl	680µl	1020µl	1.36ml
MQ Water	1.36ml	2.72ml	4.08ml	5.44ml
1.5M Tris (pH 8.8)	250µl	500µl	750µl	1ml
10% SDS	20µl	40µl	60µl	80µl
10% APS	20µl	40µl	60µl	80µl
TEMED	2µl	4µl	6µl	8µl

4. Pour the resolving gel(5ml) quickly into the caster, so as to not endup polarising the gel outside the caster.
5. To avoid bubbles, pour ethanol onto it till the brim.
6. Leave for 1 Hour.
7. Pour out ethanol.
8. Make the Stacking gel(2.5ml per cast) and pour it onto the resolving gel.

9. Put in the comb and leave for 1 hour.

The cast is ready for use.

## 13 Western Blotting

1. Blotting Buffer(50ml) :- 10ml Methanol + 35ml Water + 5ml 10 x WesternBlot Buffer
2. 6 cm x 8 cm Blotting paper/PVDF Membrane/NitroCellulose Membrane
3. 6 Blotting Sheets

Do not touch the membrane.

Set the semi-dry unit to 17V for 45 Min and 2.5A as the limit.

1. Dip the PVDF membrane in methanol(for 5 to 10 mins to activate it).
2. Throw away the methanol.
3. Add the complete(50ml) Blotting Buffer in the same container and put in the 6 blotting sheets to soak for about 6 minutes.
4. Wet the base of semi-dry unit with buffer.
5. Put in 3 sheets.
6. Then put the membrane on it with a forcep.(make sure to remove the air bubbles)
7. Cutoff the unwanted parts of the gel-cast
8. Now put hte gel-cast with protien on it.
9. Put the other 3 sheets onto the gel-cast.(again make sure to remove the air bubbles)
10. turn on the semi-dry unit and leave it for 45 minutes.  
While the protien is going from gel to the membrane make the Blocking solution.  
Blocking Buffer(100ml):- 10ml 10xTBS + 90ml MQ + 3gm SM Powder + 100µl Tween-20
11. When done with 45 minutes take the membrane out and keep with Blocking solution(100ml) for 1 hr on rocking table.
12. Throw away the Blocking solution

13. Put in the Primary antibody and soak for 60 minutes.
14. Throw away Primary antibody and wash with Blocking Buffer thrice for 10 minutes each.
15. Put in the Secondary antibody and soak for 90 minutes.
16. Throw away Secondary antibody and wash with Blocking Buffer thrice for 10 minutes each.
17. Develop with ECL.
18. Image in Chemidock.

## 14 Protien Purification

The following measurements are for 500ml of cells use by your descretion.

1. Centrifuge at 4 °C at 12k RPM for 5 mins and throw away the super-latent.
2. Add Protiase Inhibitor(200µl).  
So that after the cell covers are broken the protiases are not binded to the protine to be purified.
3. 1mg/ml lysozime (To break the cell wall)
4. Add 40ml of Buffer for resuspension then mix using pippet till pellet is dissolved.
5. Incubate(with rotator) in cold room for 40 mins.
6. Sonicate for 5 minutes at 60% amplitude with 3sec on and 3sec off pulsation.
7. Centrifuge at 12k RPM for 50 mins at 4°C and COLLECT the super-latent.
8. Resin Preparation This is the procedure for making 250µl of resin.
  - (a) 500µl of (slurry) resin.
  - (b) Centrifuge at 4k RPM for 1 minute at 4°C.
  - (c) Wash with Autoclaved MQ twice. (Add 1ml of autoclaved MQ then centrifuge with same configuration and remove the superlatent)
  - (d) Wash with Buffer twice. (Add 1ml of Buffer then centrifuge with same configuration and remove the superlatent)
9. Mix the resin with protien (forcefully).
10. Leave in cold room for more than an hour with rotator.



## 15 MiniPrep

MiniPrep is for separating DNA from cell sample.

1. Take 1.5ml of cell sample from testtube into the micro centrifuge tube.
2. Centrifuge at Room temperature at 6k RPM for 2 mins.
3. Discard the superlatent.
4. Repeat the above steps till the complete cell sample has been centrifuged.
5. Take out the Mini-prep box.
6. Add 100  $\mu$  l of solution 1.
7. Vortex until pellete is dissolved.
8. Add 100  $\mu$  l of solution 2.
9. Mix the contents by inverting but DO NOT vortex
10. Add 200  $\mu$  l of solution 3.
11. Mix the contents by inverting but DO NOT vortex
12. Centrifuge at Room temperature 13.2k RPM for 10 mins.  
The Superlatent contains the required DNA.(DO NOT throw it away)
13. Transfer the superlatent with a pipette to a microspin cup seated in a 2nd receptacle tube.  
Try to remove as little as precepetate as possible.
14. Centrifuge at Room temperature at 13.2k RPM for 1 min.
15. Open the cup, remove the microspin cup, discard the liquid from the receptacle tube and put the microspin cup back into the receptacle tube.
16. Add 750  $\mu$  l of washing buffer to microspin cup, close the cap.
17. Centrifuge at Room temperature at 13.2k RPM for 1 min.
18. Open the cup, remove the microspin cup, discard the liquid from the receptacle tube and put the microspin cup back into the receptacle tube.
19. Centrifuge at Room temperature at 13.2k RPM for 2 min.

20. Transfer the microspin cup to 1.5ml eppendorf and throw away the receptacle tube.
21. Add 50  $\mu$  l of 1xTE buffer(50°C to 55°C).
22. Keep for 10 to 15 mins.
23. Centrifuge at Room temperature at 13.2k RPM for 2 min.
24. Go for Nano Drop.(Second Floor).